# BEAUTIFLY

Scrub Peeling  
**Young**  
series

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**Technical data:**

<b>size:</b>	<b>175*43*14mm</b>
<b>weight:</b>	<b>95g</b>
<b>ultrasounds:</b>	<b>2KHZ</b>
<b>material:</b>	<b>ABS</b>
<b>rated current:</b>	<b>0.65A</b>
<b>rated voltage:</b>	<b>3.7V</b>
<b>rated power:</b>	<b>1.5W</b>
<b>battery capacity:</b>	<b>500mAh</b>
<b>included:</b>	<b>device, USB cable, user manual</b>

## Safety statement

Please read the instructions carefully before use. follow these safety guidelines when using device.

## Product parts

1. Ion Sensor
2. Protection cap
3. Ion+
4. Ion-
5. EMS indicator
6. EMS ultrasonic
7. Power button/Mode key
8. USB charging port


## Instructions

Turning on, off, operating the device

1. Press and hold the power button for 2 seconds to turn the device on. The default setting is Peeling lifting mode, and the LED1 indicator is always on.
2. Short press the power button to switch between the 4 operating modes in cycles as follow: Peeling lifting mode, LED1 indicator light is on → Power Charge mode, LED2 indicator flashes quickly → Ion+ import mode, LED3 indicator light is on, → Ion- import mode, the LED4 indicator light is on.
3. Press and hold the power button for 2 seconds to shut down.
4. The buzzer will tick once when the device is turned on and the mode is switched and the buzzer will tick twice when it is turned off.

## Charging

1. When charging: LED1, LED2, LED3, LED4 indicate breathing mode.
  2. When fully charged: LED1, LED2, LED3, LED4 indicator lights off.
- The buzzer ticks twice when the battery is low, and the indicator light shuts down.
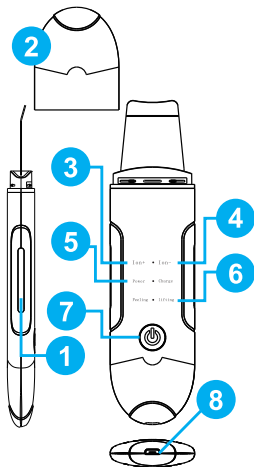
## Description of modes:

**Ion+** - positive ion export mode

Removes dirt from deep pores through positive ion export.

**Ion-** - negative ion introduction mode

Introduces active substances contained in cosmetics deep into the skin.

**Power Charge** - EMS pulse mode

Lifts and tightens, awakens sleeping collagen tissue, and provides a rejuvenating effect.

**Peeling lifting** - ultrasonic cleaning mode

High frequency ultrasonic vibration cleans the bottom of the skin, softens and removes the stratum corneum.

## Instructions for use

### Cleansing the epidermis (Ultrasounds)

1. Before starting the treatment, thoroughly cleanse the skin and remove makeup.
2. Then moisten the skin of the face with a cotton swab or an atomizer. For this purpose, you can use saline, tonic, water or hydrolate. During the treatment, the skin must be moist all the time. The treatment cannot be performed on dry skin.
3. Turn on the device and start the SONIC mode.
4. Gently move the spatula from bottom to top, lengthwise facial contours at a 45-degree angle. Repeat the action carefully, several times over the entire face.

### Cleansing the epidermis with ionization (Ultrasound + Ionization+)

1. Before starting the treatment, remember to keep the skin moist all the time. For this purpose, you can use saline, tonic, water or hydrolate. The treatment cannot be performed on dry skin.
2. Turn on the device and enter ION+ mode.
3. Gently move the spatula from bottom to top, lengthwise facial contours at a 45-degree angle. Repeat the action carefully, several times over the entire face.
4. Wipe your face with a tonic.

### The introduction of nutrients with ultrasound (Ultrasound + Ionization-)

1. After cleansing the skin of the face, apply a skin care product (cream, serum, ampoule, gel).
2. Turn on the device and start the ION mode.
3. Gently move the spatula from the bottom up along the contours of the face. The spatula should be placed in the opposite direction to that in step one. Repeat the action carefully, several times over the entire face.

If you have any doubts how to properly hold the device, please visit our channel on the YouTube platform.

 /beautifyyou


## EMS (Ultrasound + EMS currents)


1. After caring for your facial skin, turn on the device and enter the EMS mode.
2. To improve skin elasticity, move the spatula in an upward motion along the contours of your face at a 45-degree angle. Repeat the action carefully, several times over the entire face.
3. The device will turn off automatically after 10 minutes of operation
4. After the device is finished, wipe it with a cotton cloth damp cloth.

**1**  
Cheek cleaning  
care

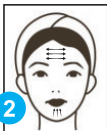**2**  
Forehead and jaw  
cleaning care

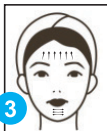**3**

**4**

Nose bridge  
cleaning care

## Tips

1. When using, hold the conductive strip on the side of the machine with your hand, so that the ion import and export function can work properly.
2. When in use, the shovel should be operated at an angle of about 45 degrees, as shown in the figure.
3. Do not stop the device in one place during use, but move it slowly.
4. If there is inflammation on the face, wait for the symptoms to disappear before using.
5. It is recommended to use 2-3 times a week. Do not use it multiple times in one day.
6. During operation, the shovel blade and the body should be clean and protected.

## Cleaning and maintenance

1. Protect the instrument from moisture, and do not place it in a dishwasher or microwave.
2. Turn off the machine before cleaning.
3. Do not use sharp, corrosive, or volatile thinner cleaning materials.
4. Use a soft cloth or cotton pad moistened with cold or warm water to wipe gently.
5. Do not rinse the device directly with water.
6. If the device is not used for a long time (about 3 months), take out and fully charge it before storing to avoid damage to the battery inside the machine.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health detriment or environmental pollution during further disposal or recycling of the equipment.

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Importer declares, on its sole responsibility, that the Scrub Peeling Young product complies with the European Parliament and Council Directive 2014/30/EU of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility, and with the European Parliament and Council Directive 2011/65/EU of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, as well as with the Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to European Parliament and Council Directive 2011/65/EU with regard to the list of restricted substances. The full text of the declaration of conformity is available at:

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