**neno®****medic T02**  
AET-R161

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## medic T02 AET-R161

USER'S MANUAL

Dear Customer,  
thank you for purchasing NENO MEDIC T02 AET-R161.  
The device you received is a non-contact infrared, reusable thermometer. The device can be used both in home environment and by professionals. The thermometer measures body temperature of adults and children that are over 3 months old.  
*Before using the device, please refer to the instructions below.*  
Keep this manual so you can refer to it later.


### 1. INTRODUCTION

- Try the device on yourself to understand how it works.
- To avoid the temperature of the surroundings influencing the measurement, the thermometer should be kept at room temperature (16-35°C) for at least half an hour before the measurement is taken.
- The person whose temperature is measured should stay inside a room with moderate temperature (16-35°C) in it for at least 20 minutes before taking the measurement.
- After doing physical exercise, wait for 30 minutes before measuring the temperature.
- When using the thermometer for a child that has its temperature measured by the thermometer and make sure that the hair does not obscure the path between the forehead and the sensor.
- If you take multiple measurements in a row, wait 5 seconds between each measurement and take them from different distances to achieve the most accurate results.
- Make sure that the sensor is clean before each measurement.
- The standard body temperature, both for adults and for children, is about 36.6 degrees Celsius. Consult your doctor if the measured temperature is much higher or lower than that.

### 2. HOW TO USE

- The thermometer measures the temperature of a person's forehead and operates at a distance of 1.5 cm to 5 cm.
- Turn on the thermometer by pressing the **I**. When you hear a beep and --°C is displayed on the screen, the unit is ready to use.
  - Place the thermometer at an appropriate distance from the forehead of the person that is being measured, with the sensor pointing towards the forehead. Press **I** and wait a second for a short beep. The temperature measurement will be shown on the display screen. If you want to perform multiple tests in a row, wait 5 seconds between measurements. This procedure can be repeated as many times as you want.
  - To turn off the device, hold **I** for about 5 seconds. The thermometer will turn off automatically after a minute of inactivity.
  - If the measurement is below 32°C, the "Lo" message will appear on the display and the thermometer will sound two beeps.
  - If the measurement is above 37.8°C, the thermometer will sound six beeps.
  - If the measurement is above 42.2°C, the "Hi" message appears on the display and the thermometer sound two beeps.

### 5. DESCRIPTION OF THE DEVICE
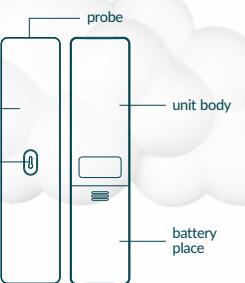
### 6. DISPLAY SCREEN


- MEASURING POSITION** – indicates the location (patient's forehead) that should be aimed at with the thermometer.
- TEMPERATURE DATA** – shows temperature measurement taken.
- TEMPERATURE UNIT** – shows whether the temperature is measured in Celsius or Fahrenheit degrees.
- BATTERY SYMBOL** – appears when batteries need to be replaced with new ones.
- MEMORY SYMBOL** – shows that the measurement currently being viewed is an archival measurement stored in the device's memory.

### 3. CHANGE THE UNIT OF MEASUREMENT


- To change the unit of measure from Celsius to Fahrenheit, hold down the **I** button for 8 seconds while the thermometer is off until the screen displays --°C. Press **I** to change the measurement unit. The thermometer will turn off automatically after 4 seconds.

### 4. MEASUREMENT MEMORY

- Make sure the thermometer is turned off.
- Hold **I** for 4 seconds to turn on memory mode. The last measurement made by the thermometer will be displayed.
- Press the **I** button to switch between measurements saved in device memory. The higher the number displayed before the temperature, the older the measurement displayed.
- Thermometer can show 32 last measurements.
- In memory mode, the thermometer turns off automatically after 12 seconds of inactivity.
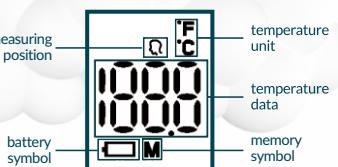
### 7. BATTERY REPLACEMENT

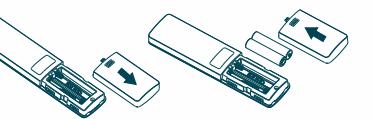- When the battery symbol appears on the display screen, immediately replace the batteries.
- Press your thumb on the lid of the battery compartment and slide it back to open the compartment between the thermometer and the device case. Pull the lid up diagonally to remove it.
- Remove the old batteries and place two new AAA batteries inside. Pay attention to the battery poles so that you don't place them the wrong way.
- Place the battery compartment cover back in place.

**Note:** Remove the batteries from the device if you will not be using it for an extended period of time. Follow local laws and regulations regarding proper battery disposal.

### 8. SAFETY MEASURES

- Do not use the thermometer in rooms where the temperature is below or above 35°C or where the humidity is below 15% RH or above 80%RH. These conditions may cause measurement errors.
- Do not expose the thermometer to temperatures below -20°C or above 65°C or humidity that is below 15%RH or above 95%RH. These conditions may damage the device.
- If the device sensor is damaged or loose, contact service immediately.
- Do not touch or blow on the sensor. This may lead to incorrect measurements.
- Do not make a self-diagnosis with the resulting measurements. If the results are of concern to you, consult a doctor.
- If the sensor is dirty, gently wipe it with soft material. Cleaning the sensor with paper towels or toilet paper can scratch and damage it.
- Incorrect installation of the battery in the device will lead to excessive battery heating and damage to the device.
- Do not expose the device to liquids. The device is not waterproof.
- The sensor is made of a special type of glass, which is difficult to damage. If you notice any damage, stop using the device immediately and contact the retailer.

### 10. MAINTENANCE, STORAGE AND CALIBRATION OF THE DEVICE

- Maintenance:**
- Remove stains from the device using soft, dry cloth.
  - Gently wipe off the stains on the device with a delicate material, use a cotton swab or cloth with alcohol to gently wipe the sensor and display. After use, the product can be sterilized with a piece of cotton or a delicate cloth soaked in medical alcohol (75% alcohol content).
  - Do not wash the appliance with water or cleaning agents containing corrosive agents or benzene.
  - Do not submerge in any liquids.
  - The device must not be stored in places where:  
thee temperature is too high or low, the moisture is too high or too low, the device is exposed directly to sunlight, near the flow of electric current or in dusty rooms. These conditions can damage the device or cause inaccurate measurements.
  - The device does not require periodic calibration. In case of measurement errors, please contact the seller.
  - The person using the thermometer should not touch the thermometer and the patient at the same time. This can lead to measurement errors.
  - Do not maintain or clean the device while it is in use.
  - The patient can measure and replace the batteries by themselves following the instructions in this manual.
  - The patient can maintain the device on their own.

### 9. ERROR MESSAGES

message	problem	solution
Hi	Temperature above 42.2°C	Use the thermometer only under the conditions specified in the manual. If the error persists, consult the service point.
Lo	Temperature below 32°C	Use the thermometer only under the conditions specified in the manual. If the error persists, consult the service point.
Err	Sensor not attached	Contact service
ErE	EEPROM memory isn't functioning properly	Contact service
	Discharged batteries	Replace the batteries with a new set
ErH	Ambient temperature too high	Go to a room with a temperature of 16-35°C
ErL	Ambient temperature too low	Go to a room with a temperature of 16-35°C

### 11. SPECIFICATION

**Measurement mode:** Forehead  
**Measurement units:** Celsius (°C) or Fahrenheit (°F)  
**Operating temperature:** 16 to 35°C  
**Operating humidity:** 15 to 80%  
**Storage temperature:** -20°C to 55°C  
**Storage humidity:** 15 to 93%  
**Measurement distance:** 15 to 50mm from the measuring point  
**Measurement range :** 32 to 42.2°C  
**Measurement accuracy :** +/- 0.3°C in 35-42°C; +/- 0.3°C in other ranges  
**Memory:** 32 last measurements  
**Dimensions:** 150x37x16.3mm  
**Weight:** 51g (without batteries)  
**Batteries:** 2xAAA (DC 3V)  
**Auto power off:** After 1 minute of inactivity  
**Product Life:** 5 years  
**Measure site:** Forehead  
**Reference site:** Armpit

### 12. EXPLANATION OF SYMBOLS

The following symbols may appear in the manual, on the thermometer, or on the rest of the set.

	Marks information that is mandatory to learn before you start using your device.
	Follow the instructions.
	BF-type part application
	Warning: Read the included documentation.
	Do not dispose of the product in a mixed municipal waste container. Dispose of the product in accordance with the disposal guidelines for electronic devices of this type.
	Transport packaging should be kept away from rain.
	Transport packaging should not be exposed to sunlight.
	Indicates the side of the transport packaging that should be facing upwards.
	The contents of the transport packaging are delicate, the packaging should be treated with extreme care.
	Indicates the temperature range at which the transport packaging should be stored.
	Non-ionizing electromagnetic radiation
	Don't roll.
	LOT Number
	Production date
	The product cannot be used after the date listed next to this icon.
	The product is resistant to solid objects with a diameter of 12.5mm or more and against vertically falling water droplets when the device is raised at an angle of up to 15 degrees.
	Manufacturer
	CE mark: The product complies with the Medical Device Directive 93/42/EEC
	Authorized representative in the European Community.

### 13. ELECTROMAGNETIC COMPATIBILITY

The ME DEVICE is adapted for use in home and similar environments for healthcare purposes.

**Warning:** Do not start the device near high-frequency surgical devices or in magnetic resonance rooms with RF shielding, as there are high amounts of electromagnetic disturbances in these locations.

**Note:** Using the device near other equipment may lead to measurement irregularities and should be avoided. If the use of two devices in each other's vicinity is unavoidable, both devices should be carefully monitored to ensure that they are working properly.

If any: a list of all cables and maximum cable length (if applicable), transducers and other ACCESSORIES which are replaceable by THE RESPONSIBLE ORGANIZATION and which may affect the compliance with the standards of ME EQUIPMENT or ME SYSTEM with the requirements set out in Clause 7 (EMISSIONS) and Clause 8 (IMMUNITY). ACCESSORIES can be specified in general (e.g. shielded cable) or specifically (e.g. specifying the MANUFACTURER and TYPE OF EQUIPMENT).

If any: the functioning of medical equipment or medical system has been defined as ESSENTIAL PERFORMANCE and a description of the changes that the OPERATOR may experience if ESSENTIAL PERFORMANCE is disturbed by EM DISTURBANCES (the term ESSENTIAL PERFORMANCE does not need to be used).

1. All necessary instructions for maintaining BASIC SAFETY and ESSENTIAL PERFORMANCE with regard for EM disturbances for the period of the device's expected service life.

2. Guidance and manufacturer's declaration - electromagnetic emissions and immunity.

Table 2

Guidance and manufacturer's declaration - electromagnetic immunity

Immunity test	IEC 60601-1-2 Test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	+/- 8 KV contact +/- 8 KV air	+/- 8 KV contact +/- 8 KV air
Electrical fast transient/burst IEC 61000-4-4	Power supply cables: +/- 2KV Input/output cables: +/- 1KV	N/A
Surges IEC 61000-4-5	Line(s) to line(s): +/- 1KV Line(s) to earth: +/- 2KV Repetition frequency 100 kHz	N/A
Voltage dips, short interruptions and voltage variations in power supply cables IEC 61000-4-11	0% 0.5cycle -10%, -45%, 90%, 135%, 180°, 225°, 275° and 315° 0% 1 cycle And 70% 25/30 cycles Single phase: at 0% 300 cycles	N/A
Magnetic field frequency IEC 61000-4-8	30 A/m 50Hz/60Hz	30 A/m 50Hz/60Hz
Conducted radio frequencies IEC 61000-4-6	150KHz to 80MHz: 3Vrms 6Vrms (in ISM and amateur radio bands) 80% Am at 1kHz	N/A
Radiated radio waves IEC 61000-4-3	10 V/m 80MHz - 2.7GHz 80% Am at 1kHz	10 V/m 80MHz - 2.7GHz 80% Am at 1kHz

NOTE: Ut is the AC mains voltage before applying the test level.

Table 1

Guidance and manufacturer's declaration - electromagnetic emission

Emission test	Standard met
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class B
Harmonic emissions IEC 61000-3-2	N/A
Voltage fluctuations/flicker emissions IEC 61000-3-3	N/A

Table 3

Guidance and manufacturer's declaration - electromagnetic immunity

Test frequency (MHz)	Band (MHz)	Service	Modulation	Modulation (mW)	Distance (m)	Measured level (mV/m)
385	380-390	TETRA 400	Pulse modulation 18Hz	1.8	0.3	27
450	380-390	GMRS 460, FRS 460	FM +/- 5 kHz, 1 kHz sine	2	0.3	28
710	704-787	LTE band 13, 17	Pulse modulation 217Hz	0.2	0.3	9
810	800-960	GSM 800/900, TETRA 800, IDEN 820, CDMA 850, LTE Band 5	Pulse modulation 18Hz	2	0.3	28
1720	1700-1970	GSM 1800, CDMA 1900; DECT, LTE Band 1, 3, 4, 25, UMTS	Pulse modulation 217Hz	2	0.3	28
2450	2400-2570	Bluetooth, WLAN 802.11 b/g/n, 2.4GHz, 2.5GHz, LTE Pasmo 7	Pulse modulation 217Hz	2	0.3	28
5240	5100-5785	WLAN 802.11 a/n	Pulse modulation 217Hz	0.2	0.3	9

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