S·E·Y·S·S·O

PLANET

MERCURX/NEPTUNX
SONIC TOOTHBRUSH

PLANET

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EN / PL

Please read the user manual before use. Keep the manual for possible future reference.

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1 Safety instructions

Warning:


1. Before connecting the product to the power supply, Make sure that your outlet voltage corresponds to the voltage stated on the rating label – 100-240V~,50-60Hz.
2. Do not charge the toothbrush in damp places.
3. Do not pull the cord when disconnecting the charger. Do not touch the cord with wet hands.
4. If the charger falls into water during charging, disconnect the power supply before taking it out.
5. Do not immerse the toothbrush in water at a temperature above 40°.
6. Do not charge, use or leave the device if it is exposed to high temperatures or environmental factors.
7. The device is not intended for unattended use by children or persons with reduced motor, mental or cognitive abilities. Keep the device out of their reach. The product is not a toy.
8. Children or persons with reduced motor, mental or cognitive abilities may only use the toothbrush under appropriate supervision.
9. Do not press the toothbrush against your teeth or gums, or brush in the same place for a long time, as this can damage your gums.
10. Inspect the device thoroughly before use. Do not use it if you detect any visible damage.
11. If you have undergone any oral surgery or treatment of your teeth and gums in the last two months, ask your dentist if using the toothbrush is advisable.

12. The device may cause slight bleeding of the gums during the first days of use. If the bleeding does not stop within a week or is unusually heavy, contact your doctor.
13. Persons with an implanted cardiac pacemaker or other electronic device should consult a doctor before using the toothbrush.
14. Do not attempt to repair, modify or disassemble parts of the device by yourself. In case of damage, please contact the customer service centre.
15. We recommend replacing the brush head every 3 months or more often, depending on the frequency of use.
16. If you use toothpastes containing brightening or whitening substances, e.g. sodium bicarbonate or other bicarbonates, wash the brush thoroughly after each use.
17. The device is designed for brushing teeth, gums and tongue. Do not use it for any other purpose. If you feel any pain or discomfort after using the toothbrush, or if you have any concerns about the impact of the device on your health, contact your doctor. The power cord should not be replaced. If the cord is damaged the appliance should not be used.

Electromotive force:

The SEYSSO sonic toothbrush meets all standards for the operation of electromagnetic devices. If used properly and in accordance with this manual, it will not pose any danger to the user.

2 Toothbrush design

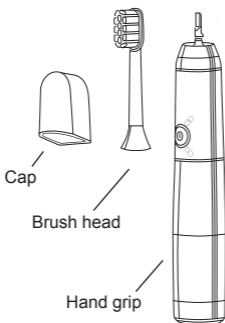3 Charging before the first use:


1. Discharge the toothbrush completely before the first charging. Next, charge the toothbrush until it is fully charged. The first charging should last 32 hours. During this time, you can remove the toothbrush from the charger to clean your teeth.
2. A LED ring on the toothbrush flashes in four different colors during charging. The LED ring goes out 2 minutes after the toothbrush is fully charged.
3. Each subsequent recharging should last 32 hours or until the LED ring goes out.

ATTENTION: Switch off the toothbrush before the first charging. Place the charging base on a flat and stable surface and align the slot on the bottom of the toothbrush with the locking pin of the charging base. Place the toothbrush on the charging base and connect the plug to a power outlet.

ATTENTION: Before charging, make sure that the output current specification is compatible with 100-240 V AC 50 Hz/60 Hz 1.4 W. Otherwise, do not connect the charger to the power supply.

4 Brushing:

Step 1: Attach the replaceable head to the hand grip.

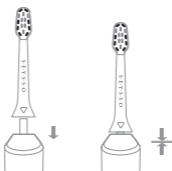* The diagram shows the correct way to attach the head to the toothbrush

Step 2: Next, apply an appropriate amount of toothpaste to the bristles of the toothbrush. We recommend using toothpaste dedicated to sonic brushes (Sonic Friendly). We recommend the SEYSSO toothpastes which, thanks to their excellent ingredients and texture, fully utilise the potential of sonic toothbrushes.

Step 3: Put the head into your mouth before switching on the toothbrush. Position the toothbrush near your teeth so that its bristles are at a 45-degree angle to your gum line and switch on the toothbrush by pressing and holding the function button for 3 seconds.

Step 4: Briefly press the function button while the toothbrush is running to select the appropriate brushing mode. Each time you press the button, the device switches between the four programs in sequence and changes the colour of the LED ring.

Mode	LED ring color	Symbol on the toothbrush	Function
Clean	Yellow	C	Teeth cleaning
White	Blue	W	Teeth whitening
Polish	Red	P	Polishing front teeth
Massage	Green	M	Gum massage and teeth cleaning

Step 5: Slowly move the toothbrush along the gum line without pressing your teeth too hard. The toothbrush has a built-in automatic 2-minute timer to help you keep the recommended brushing time. Thanks to 30-second intervals, signalled by a short pause of the toothbrush operation, you can evenly distribute the brushing time to all quarters of your teeth.

Step 6: The device will automatically switch off after 2 minutes. Press and hold the button for 3 seconds to turn off the toothbrush while it is running.

Step 7: After brushing your teeth, remove the head from the hand grip, rinse it with water and dry it before storing it. Remember to remove the head after each use of the toothbrush. When the battery runs out, place the toothbrush on the charging base and charge it according to the above instructions.

Step 8: When the battery level drops below 3.3 V, the LED ring will start to flash slowly, whereas when it drops below 3 V, the LED ring will be flashing quickly, and the toothbrush will turn off automatically.

5 Replaceable heads:

The set includes the SEYSSO Planet XCleaning brush head. Use only original SEYSSO heads dedicated to the toothbrush. Any damage caused by using non-original heads will not be covered by the manufacturer's warranty.

6 Cleaning and maintenance

WARNING Disconnect the brush from the power supply before cleaning. Do not immerse the charging base in water. Rinse the hand grip and the toothbrush head with water at a temperature not exceeding 40°. Wipe the charging base and hand grip with a soft cloth. Never use strong cleaning agents.

WARNING: Do not use chemical cleaners or tools that may damage the toothbrush. Do not wash the head, charger or hand grip in a dishwasher.

WARNING: Li-ion batteries must be charged at room temperature. Too low (below 0°) or too high (above 45°) temperature may cause serious safety problems.

7 Warranty and servicing

The SEYSSO toothbrush is covered by a two-year manufacturer's warranty, which is valid from the date of purchase. During the warranty period, the manufacturer undertakes to repair or replace a defective toothbrush with a new one if the defect is caused by a structural flaw. The warranty does not cover damage resulting from improper use of the device.

8 Quality guarantee – restrictive clause

The quality guarantee covers only SEYSSO devices.

Explanation:

Please fill in the warranty card thoroughly and keep it for possible future use. The user is entitled to a two-year period of free service, which starts on the day of purchase.

The user shall not be entitled to the free service if:

1. The product is damaged due to improper use or storage;
2. The product has been dismantled or repaired by a repair facility other than that indicated by the manufacturer;
3. The user has no valid proof of purchase of the product;
4. Product damage is caused by an unstable voltage in the mains or a voltage exceeding the specified range;
5. The warranty period for the product has expired.

SEYSSO Planet NeptunX/MercurX:
100-240 V~, 50/60 Hz; 1,4W;

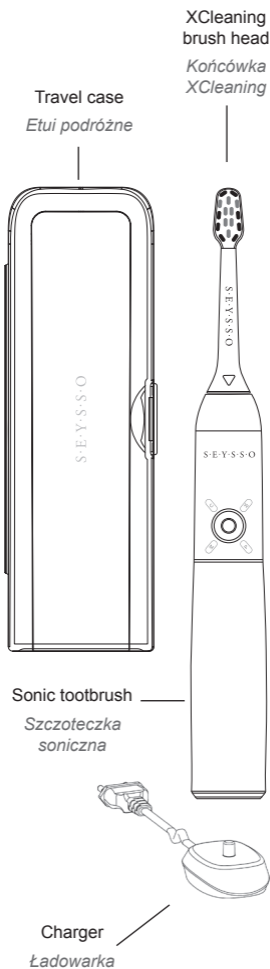
**Enviromentally
friendly disposal**

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